


Alexis™ ViU™

Wound Protector-Retractor

SETUP GUIDE



1 Plug the nonsterile, Alexis ViU reusable power supply* into a properly grounded receptacle (i.e., wall outlet or surgical boom).

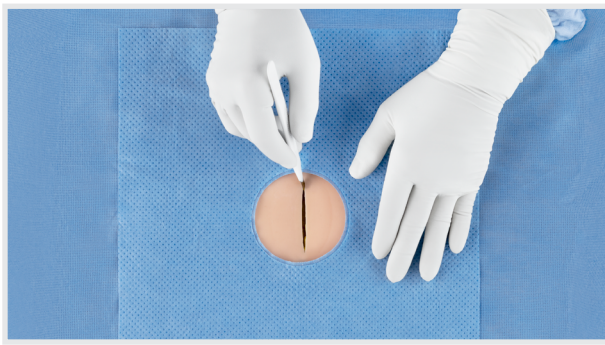

2 Make the incision at the desired location.



3 Insert the clear inner ring through the incision, keeping the power cord and tether (if applicable) outside the sheath and patient.




4 Sweep the area to ensure that no tissue entrapment has occurred.



5 Grasp the white outer ring at the 10 and 2 o'clock positions, and pull up until the clear inner ring sits against the tissue layer.



6 Roll the white outer ring inward or outward until the desired retraction is achieved. Sweep the area again to ensure that no tissue entrapment has occurred.

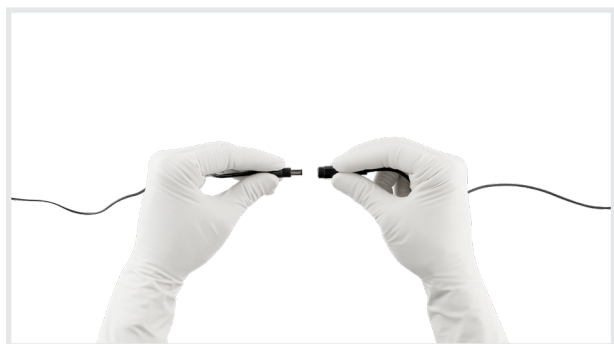

7 Connect the device power cord to the reusable power supply in the **nonsterile field**.



8 Perform the surgery through the illuminated incision site.



9 Disconnect the device from the power supply, and gently remove it by pulling the clear inner ring or tether tag (if applicable) from the incision.

ADDITIONAL INSTRUCTIONS

- The device should only be used with the Alexis ViU reusable power supply (CLEXX).†
- After use, thoroughly clean the surfaces of the reusable power supply following facility protocol.
- Store the reusable power supply according to your facility practices.

Visit appliedmedical.com/alexisviu for more information.

* The Alexis ViU reusable power supply is sold separately.

† The model number for the Alexis ViU reusable power supply varies by region to comply with local electrical standards.

The devices listed may not be approved in all markets.

Devices listed must be administered by a healthcare professional. Please contact your Applied Medical representative for more information on availability. This information is intended for dissemination exclusively to healthcare professionals and is not intended to replace labeling and Instructions for Use (IFU). Please refer to the IFU for the indications, contraindications, warnings, precautions, instructions and other information.

© 2026 Applied Medical Resources Corporation, Rancho Santa Margarita, California. All rights reserved. Applied Medical, the Applied Medical logo design, and marks designated with the symbol ™ are trademarks of Applied Medical Resources Corporation, registered in one or more of the following countries: Australia, Canada, Japan, South Korea, the United States, the United Kingdom, and/or the European Union. ViU is a trademark of Applied Medical Resources Corporation, pending registration.

03/2026 732615-EN-A