



Silhouette[®]


Ureteral Stents

Silhouette

Ureteral Stents

Revolutionary technology, innovative design and proprietary material make the Silhouette ureteral stents different from standard double J stents.

Offering all of the clinical advantages with maximum patient tolerance, Silhouette ureteral stents are the next generation in ureteral patency.

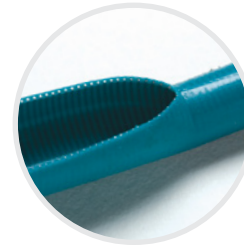


Advanced Technology

Coil reinforced wall maximizes inner lumen allowing enhanced drainage comparable to large stents, while providing unparalleled patient tolerance and effortless placement

Superior Compression Resistance

Coil-reinforced body optimizes patency by minimizing kinking and buckling and maximizing compression resistance ^{1,2,3,4}

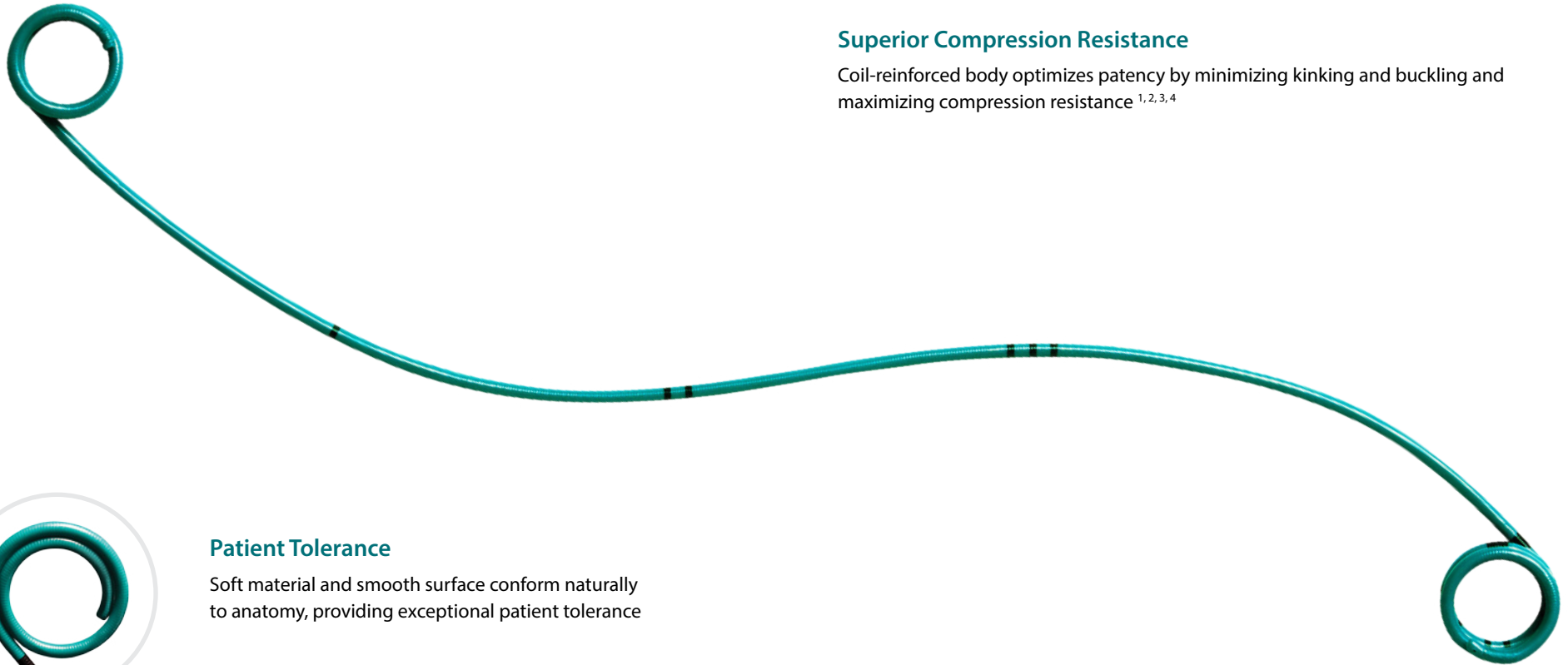

Patient Tolerance


Soft material and smooth surface conform naturally to anatomy, providing exceptional patient tolerance

Easy Insertion

Hydrophilic coating and superior column strength provide unprecedented ease of insertion and advancement

MRI Compatible

Silhouette ureteral patency devices have been tested and approved as MR-conditional (3-Tesla or less) for use in magnetic resonance imaging

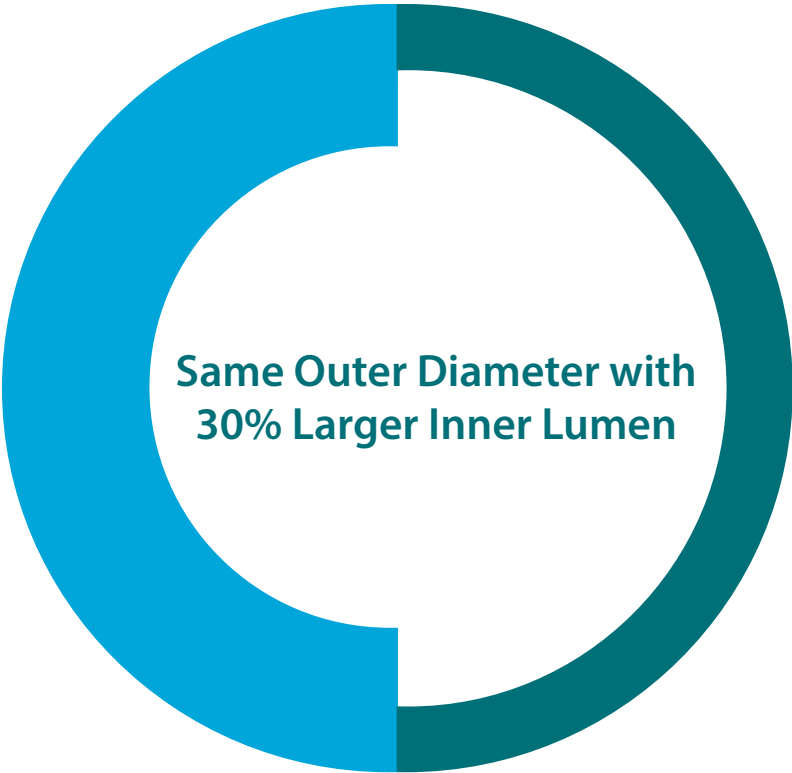

Standard 6F
Stents all lengths

Silhouette XtraFlo 6F
Stents all lengths



Silhouette Pediatric ureteral stents

Reorder No.	Size	Length	Qty/Box
B3883*	4F	8cm	1/Box
B3884*	4F	10cm	1/Box
B3885*	4F	12cm	1/Box

* Includes .035" UroWIRE® guidewire



Silhouette Comfort ureteral stents

Reorder No.	Size	Length	Qty/Box
B3829	4.6F	10cm	1/Box
B3830	4.6F	12cm	1/Box
B3831	4.6F	14cm	1/Box
B3832	4.6F	16cm	1/Box
B3833	4.6F	18cm	1/Box
B3834	4.6F	20cm	1/Box
B3835	4.6F	22cm	1/Box
B3836	4.6F	24cm	1/Box
B3837	4.6F	26cm	1/Box
B3838	4.6F	28cm	1/Box
B3839	4.6F	30cm	1/Box
B3827	4.6F	22-30cm	1/Box



Silhouette XtraFlo ureteral stents

Reorder No.	Size	Length	Qty/Box
B3849	6F	10cm	1/Box
B3850	6F	12cm	1/Box
B3851	6F	14cm	1/Box
B3852	6F	16cm	1/Box
B3853	6F	18cm	1/Box
B3854	6F	20cm	1/Box
B3855	6F	22cm	1/Box
B3856	6F	24cm	1/Box
B3857	6F	26cm	1/Box
B3858	6F	28cm	1/Box
B3859	6F	30cm	1/Box
B3847	6F	22-30cm	1/Box



Silhouette Scaffold ureteral stents

Reorder No.	Size	Length	Qty/Box
B3870	8F	14cm	1/Box
B3871	8F	16cm	1/Box
B3872	8F	18cm	1/Box
B3873	8F	20cm	1/Box
B3874	8F	22cm	1/Box
B3875	8F	24cm	1/Box
B3876	8F	26cm	1/Box
B3877	8F	28cm	1/Box
B3878	8F	30cm	1/Box
B3869	8F	22-30cm	1/Box

1. Christman M., MD., et al: Analysis of Ureteral Stent Compression Force and its Role in Malignant Obstruction. *J Urol.* 2009 Jan; 181 (1): 392-6. Epub 2008 Nov 17.

2. Miyaoka R., MD., et al: Resistance to Extrinsic Compression and Maintenance of Intraluminal Flow in Coil-Reinforced Stents (Silhouette Scaffold Device): an in Vitro Study. *J Endourol.* 2010 Apr; 24(4):595-8.

3. Christman M., MD., et al: Analysis of Ureteric Stent Kinking Forces: The Role of Curvature in Stent Failure. *BJU Int.* 2010 Mar; 105(6):866-9; discussion 868-9. Epub 2009 Aug 27.

4. Pedro, R., et al: Wire-Based Ureteral Stents: Impact on Tensile Strength and Compression. *Urology* 70(6) 2007; 1057-1059.

Visit www.appliedmedical.com/silhouette
for more information

© 2017 Applied Medical Resources Corporation. All rights reserved.

Applied Medical, the Applied Medical logo design and marks designated with a © are trademarks of Applied Medical Resources Corporation, registered in one or more of the following countries: Australia, Canada, Japan, South Korea, the United States and/or the European Union. SC01387J

